


Healthy Massage Chair
ISO9001:2000

imedics
www.imedics.com.au

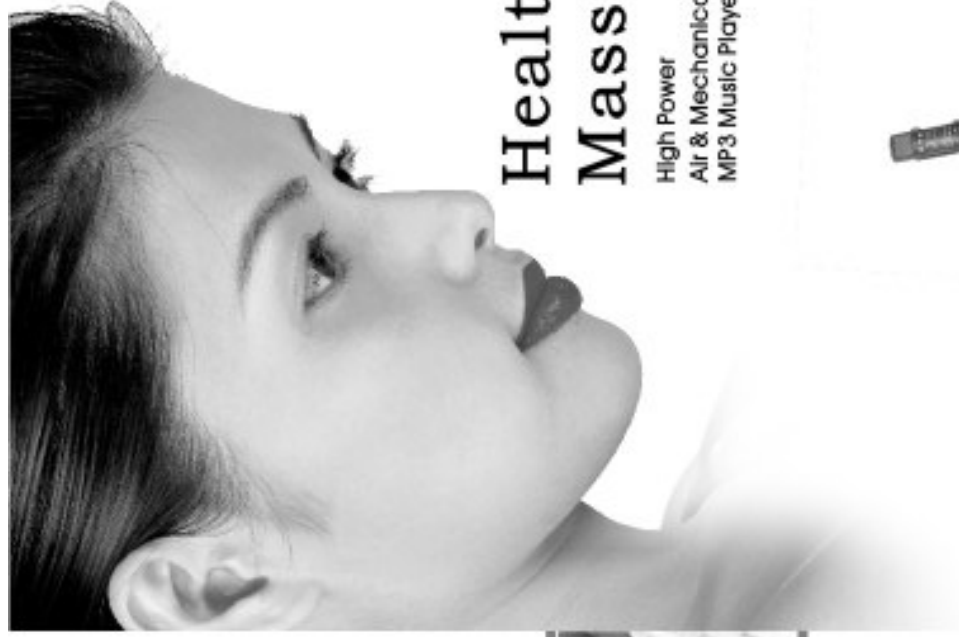

Healthy Massage Chair

High Power
Air & Mechanical Massage Combination
MP3 Music Player



ALL RIGHTS RESERVED

 You can just email us your
warranty card online
warranty eWarranty card
www.imedics.com.au



MODEL: AMBIENCE

USERS MANUAL

TYPE: AMBIENCE

⚠ Caution:



- Before using the chair, lift the backrest pad to check for tears on the backrest centre cloth and its surrounding upholstered areas.
- Do not use the chair if the backrest centre cloth or its surrounding upholstered areas of the chair is torn, regardless of the severity of the tear(s).
- Remove the AC plug from the electrical outlet immediately and contact Medica Service centre. Continuous usage of the chair in this condition may cause injury to the user or lead to an electric shock.

TYPE: AMBIENCE

⚠ Caution:

- The following users must consult their doctor before use:
- Those who are under medical care or are feeling unwell.
- Those who have spinal disorders, bent or abnormal spinal conditions caused by past accident or illness.
- Those who would like to apply the massage to an area that was or is being medically treated.
- Those who are suffering from malignant tumours.
- Those who are suffering from heart diseases or severe diabetes.
- Those who are pregnant or menstruating.
- Those who are under medical rest, or rest-cure ordered by doctors.
- Those who have osteoporosis.

- Do not operate the chair when the AC plug is loose, the power cord is damaged, or the electrical outlet is excessively loose.
- Do not wind the AC cord or the remote control cord around the chair, or place them under the chair. Do not damage, modify or handle the AC cord, Headphone cord or the remote control cord roughly.
- Do not allow children or pets to play around the chair, i.e. behind the backrest, under the seat or footrest of the chair.
- Ensure that there is no obstruction in the paths of the backrest and the footrest during the adjustment of the recline angle of the chair.
- Do not use the chair if your body is wet. Do not operate the chair with wet fingers.
- Do not stand on the chair. Do not sit on any other part of the chair except the seat, to prevent damage to the chair.
- Please read the operation manual thoroughly for the correct methods of usage, especially the notes on safety.

SPECIFICATIONS

Model: AMBIENCE
Voltage: 240V ~
Frequency: 50Hz
Power: 220W
Safe Design: Class I

Before use, please carefully read this users manual .Pay attention to the notice terms in this manual.


After carefully reading, please keep this users manual within reach at all times for random help.

AFTER USE.....

- Keep the remote control in the bracket
- To avoid any accidents or damage, make sure to turn off power and remove plug from power socket.

Storing way:

Firstly reposit footrest and backrest to the lowest and upmost position separately, turn off all massage functions and remove plug from power socket, then release lock lever as figured to lay the backrest onto the seat cushion.

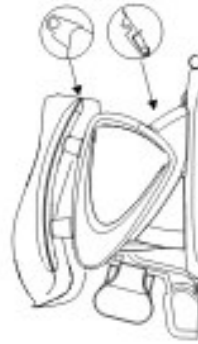

Deluxe Massage Chair IMPORTANT SAFETY INSTRUCTIONS

TERMS FOR USING

	When descending the footrest or backrest, use moderate pressure and be sure there are no obstructions or small children near or under the footrest.
	Before use, uncover back cushion and make sure inside cushion is in good order as with all other parts. Should any part be damaged contact distributor immediately and do not use until repaired.
	To avoid any fire or damage to the product, rated voltage must be used only.
	If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or similarly qualified person in order to avoid hazard.
	Do not allow children or unconscious people to use product without any assistance.
	Do not allow children to play on the product.

	Make sure to load not more than two people in case any hazard or injury.
	Do not use pointed or sharp objects on or near product at any time and avoid dragging on material of product.
	Do not use it while smoking to avoid fire.
	Do not use other medical apparatus when operating the product.
	Should any abnormalities occur when using, stop to consult with doctor.
	When disconnecting plug from power socket, be sure to hold insulating bush instead of power cord.
	Do not use petrol, caustic agents and spray pesticide for cleaning.
	Do not sit on backrest, armrest and footrest to avoid any injury or damage.
	Allow only authorized technicians to handle any breakdown, repair or alteration to the product.
	When staining or maintenance, remove plug from power socket and make sure not to operate by wet hands to avoid any hazard.
	After use, turn power off.


OPERATING METHOD




1. The stimulated backrest mechanism design to the back shape.

2. Delicate massage presents humanistic and overall massage enjoyment.


The moving backrest massage mechanism between neck and waist performs kneading, tapping, pressing and synchronization massage modes. Meanwhile adjustable speed and width help to relieve muscular fatigue and increase blood circulation.

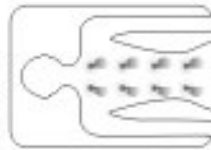


Knead

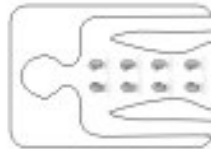

Tapping

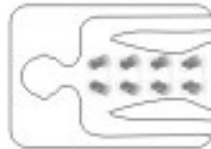


Synchronization

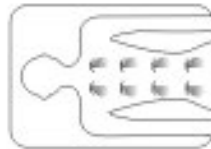

Pressing

3. Lower back adopts wavelet whole air pressure massage

Air bags at thigh, buttock and leg parts have been inflating and deflating air alternatively to offer repeated air pressure massage and really helps to relax nerve and improve blood circulation for lower body.

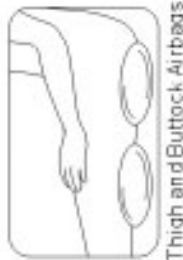

Thigh and Buttock Airbags



Legs and Feet Airbags

4. Overall Automatic Programs

- Stop all massage functions and reposition backrest and footrest automatically
- A** Upper body AUTO: Select upper body (the whole backrest) automatic programs
- B** Neck & Shoulder: Select neck and shoulder parts automatic programs
- C** Back & Waist: Select back and waist parts automatic programs
- D** Airbags AUTO: Turn on airbags automatic programs
- E** Air Pressure Strength: Select air pressure strength from "weak" and "strong"

MP3 MUSIC PLAYER

The product adopts implanted MP3 music player, either speaker or headphone are available to enjoy music

- Music renewal

① Please make sure to shut off power before following operations

- 1 Plug the memory stick (store songs) in the USB groove of the computer, with install windows 2000 or windows XP.
- 2 Computer will automatically recognize MP3 music player and install device driver.
- 3 There will appear a driver letter.
- 4 You can then copy, delete or paste etc. operations on this disk.
- 5 After renewal, please be sure to carry out "safely delete" operation.
- 6 Plug in the USB jack of the massage chair, you can enjoy the renewed music.

For further instructions, please refer to "music play control" term from "NAME AND FUNCTION FOR EACH PART"

STOP AND STORING

- Setting time is up or press "power" button, backrest, footrest all reposition to original positions.



- Press "settle" button, not only carry the same functions as "power" button, but backrest, footrest all reposition to original positions.



POWER ON

- Connect plug to power socket
- Switch on (located at lower back cover)

TIME SETTING


The setting time is massage time, it shows on LCD by surplus time. When setting time is up, all massage functions stop, and massage mechanism reposition. Setting time can be selected from 10, 20 or 30 minutes and default time is 20 minutes.

BACKWARD RECLINING

Keep pressing "▲" button, backrest gradually goes down while footrest comes up. When it reaches desired angle, release "▲" button, backrest and footrest stop moving. When backrest and footrest reach upmost and lowest position, a reminding sound occurs.




Backward reclining

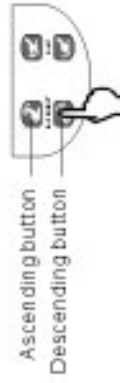


FORWARD RECLINING

The operating method is the same than "backward reclining".

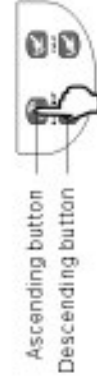


Forward reclining



External Structure

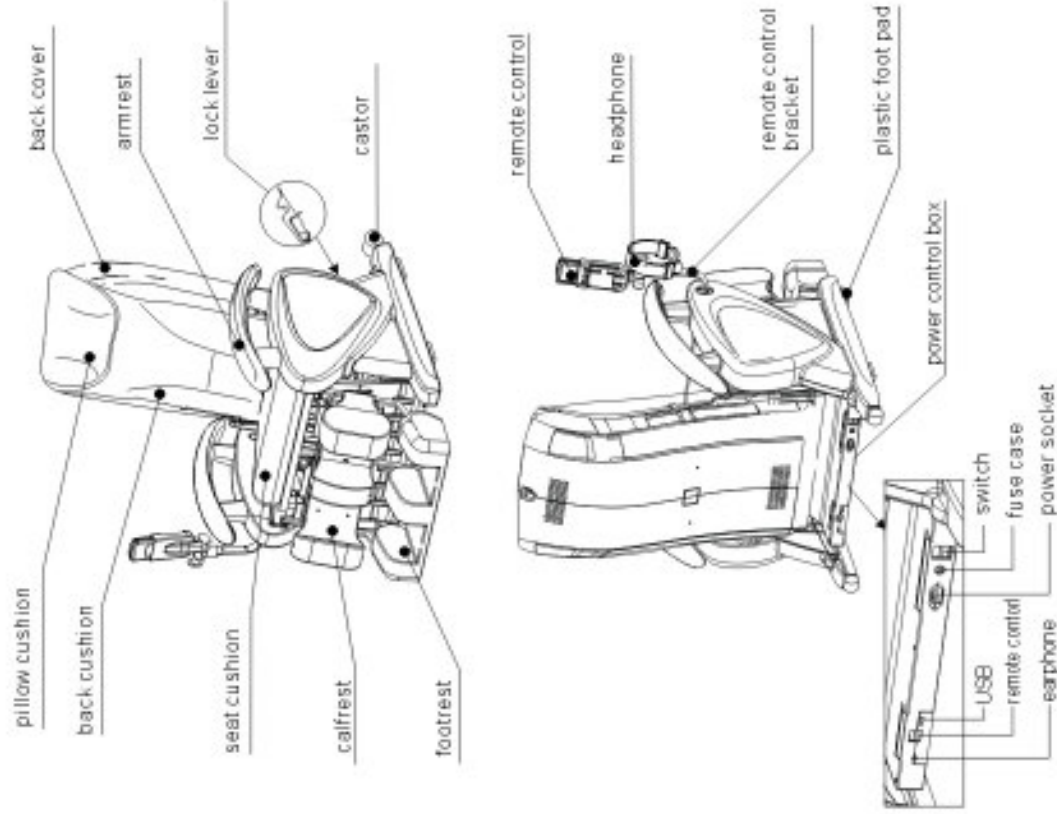

INTERNAL STRUCTURE

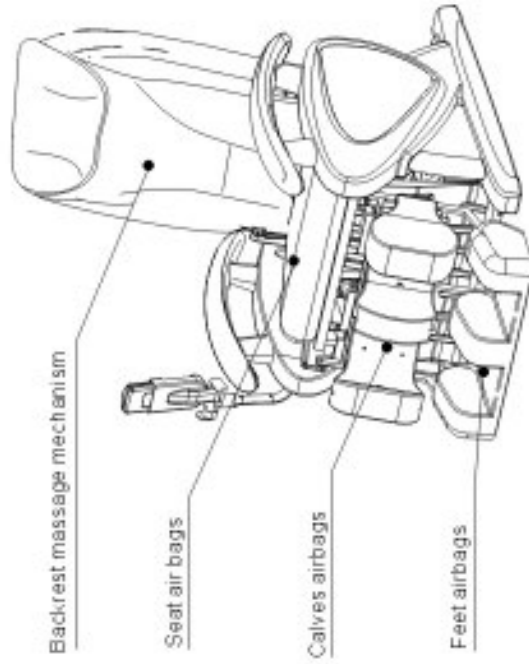

INSTALLATION METHOD

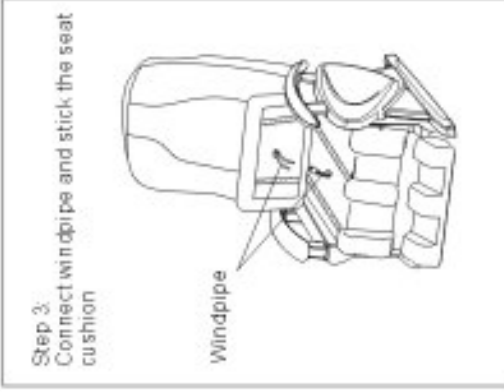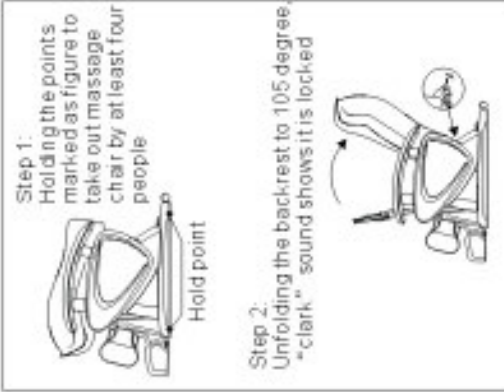


PREPARATIONS BEFORE USE

INSTALLATION PLACE

Safety distance around massage chair.

Horizontal interval of backrest between upmost and lowest point should be more than 50cm.

Horizontal interval of footrest between upmost and lowest point should be more than 20cm

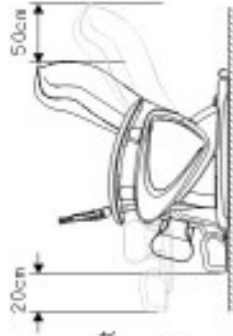

FLOOR PROTECTION

The product could damage the floor, so a medium material under the chair (carpet for instance) is strongly advised to use for floor protection



SURROUNDINGS



Don't keep the product in a bathroom or damp place to avoid electricity leaking

Don't keep the product near stove or heating devices to avoid fire or damage to the leather cover

GROUNDING CONNECTION



The product belongs to classification I electric appliance. Please make sure to confirm ground connection to avoid electricity leaking or any accident when using


MOVEMENT WAYS

1. Lift the front of the chair by using footrest and push with two hands to start moving.
2. When moving chair do not have persons or objects on chair to avoid overturning.
3. Do not move or drag the massage chair on a delicate floor.
4. Do not move chair by holding on to the seat or backrest, footrest must be used only.
5. The chair is heavy, so handle with care.



NAME AND FUNCTION FOR EACH PART

OUTER CONTROL PANEL

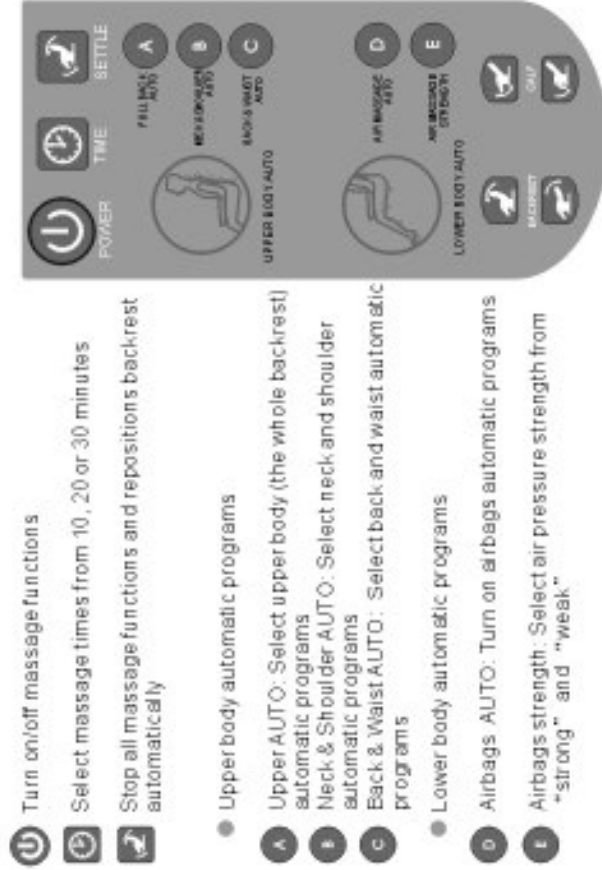

- Turn on/off massage functions
- Select massage times from 10, 20 or 30 minutes
- Stop all massage functions and reposition's backrest automatically

- Upper body automatic programs
 - Upper AUTO: Select upper body (the whole backrest) automatic programs
 - Neck & Shoulder AUTO: Select neck and shoulder automatic programs
 - Back & Waist AUTO: Select back and waist automatic programs
- Lower body automatic programs
- Airbags AUTO: Turn on airbags automatic programs
- Airbags strength: Select air pressure strength from "strong" and "weak"

- Backrest ascending: When backrest comes up, footrest will go down in the meantime. Release button, it stops moving. When backrest reaches upmost, a reminding sound occurs.

- Backrest descending: Keep pressing the button, backrest goes down and footrest comes up in the meantime, release it, it stops moving. When backrest goes down to the lowest, a reminding sound occurs.


- Footrest ascending: Keep pressing the button, footrest comes up, release it, footrest stops ascending. When footrest reaches upmost, a reminding sound occurs.

- Footrest descending: Keep pressing the button, footrest goes down, release it, footrest stops descending. When footrest reaches to the lowest, a reminding sound occurs.

INNER CONTROL PANEL

music play control

- ▶ Play : play music
- ▶ Pause: music playing pause
- ▶ Up: select previous song
- ▶ Down: select next song
- ▶ Turn up volume
- ▶ Turn down volume

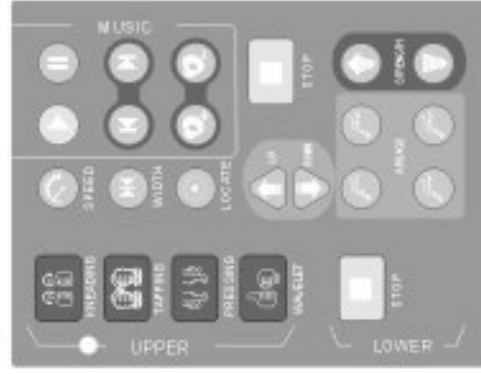


Lower Body Massage

- Airbags massage manual control
- Select soles and calves airbags massage
- Select soles and seat airbags massage
- Select calves and seat airbag massage
- Select soles, calves and seat airbags massage
- Air pressure strength adjustment (valid for both AUTOMANUAL modes)
- Select "strong" air pressure massage
- Select "weak" air pressure massage
- Turn off air pressure massage (valid for both AUTOMANUAL modes)
- Turn off air pressure massage (valid for both AUTOMANUAL modes)


Upper body massage

- ▶ Knead: select kneading massage mode
- ▶ Tapping: select tapping massage mode
- ▶ Pressing: select shiatsu massage mode
- ▶ Synchronization: select knead/tapping Synchronal massage mode
- ▶ Speed adjustment: Select three massage speeds from SLOW, MODERATE or FAST. This adjustment is invalid to "Pressing" massage mode.
- ▶ Width adjustment: Select massage width from NARROW, MIDDLE or WIDE. This adjustment is only valid to "Tapping" and "Pressing" massage modes
- ▶ Locate: Select massage positions from OVERALL, PARTIAL (shoulder, backrest and waist) or POINT. OVERALL refers to the whole back; PARTIAL refers to a partial massage area; POINT refers to massage a certain point, under this mode, micro-adjustment is optional.
- ▶ Under POINT mode, select point's upwards
- ▶ Under POINT mode, select point's downwards
- ▶ Stop bedrest massage functions and reposition back massage mechanism


LCD

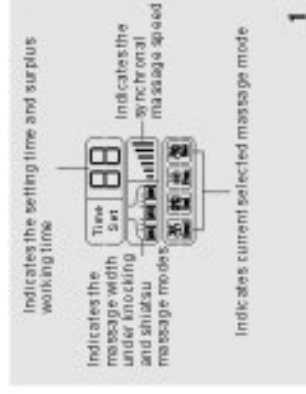

- Indicates the setting time and surplus working time
- Indicates the message width under knocking and shiatsu massage modes
- Indicates current selected massage mode

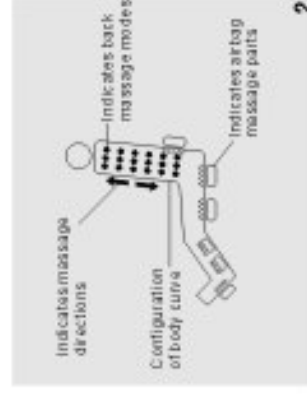


- Indicates massage directions
- Indicates back massage modes
- Indicates airbag massage parts

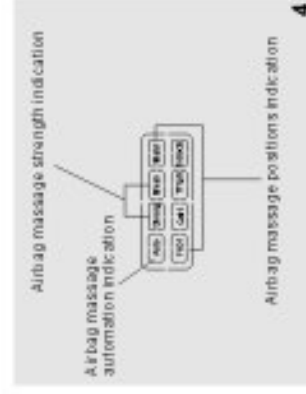

- Airbag massage strength indication
- Airbag massage automation indication
- Airbag massage parts indication



- Indicates back massage positions
- Indicates back massage modes

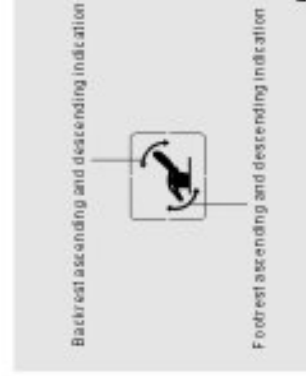

- Backrest ascending and descending indication
- Footrest ascending and descending indication